

Class II Type B1 Laminar Flow Biological Safety Cabinet Ordering Criteria for the National Institutes of Health

Background

This document provides guidelines for the design and installation of Class II, Type B1 Biological Safety Cabinets (BSCs) at National Institutes of Health facilities.

Biological Safety Cabinets

A Class II Biological Safety Cabinet (BSC) is commonly found in laboratories working with infectious agents or with tissue culture. By definition, a Class II BSC is a ventilated cabinet, which provides personnel, product and environmental protection. It has an open front with inward airflow for personnel protection downward High Efficiency Particulate Air (HEPA) filtered laminar airflow over the work surface for product protection and HEPA filtered exhausted air for environmental protection. In order for the BSC to operate as designed and be properly maintained, the following guidelines have been developed.

Before purchasing any BSCs for installation in an NIH facility, contact the Division of Occupational Health and Safety (DOHS) for plan reviews, work reviews, and approval at 301 496-3457.

Locating the Biological Safety Cabinet

- Locate the BSC out of the traffic pattern and away from room air currents that could disrupt the air barrier at the work access opening.
- Do not locate the BSC where room supply air inlets blow across the front opening or onto the cabinet's exhaust filter.
- Provide a minimum of 4 inches clearance on the rear and the non-utility side, and 6 inches clearance on the utility side of all cabinets and locate the power receptacle high on the wall so that the unit may easily be unplugged for servicing.
- To help assure that the 6-inch clearance provides enough space for utility connection, use flexible gas connections.
- Seismically secure all BSCs using a connection that may easily be disconnected for servicing and relocation.
- See Drawing – Class II Type B1 Biosafety Cabinets (page 4).

Venting Class II Type B1 Biological Safety Cabinets:

- Direct-duct all Class II, Type B1 BSCs to an exhaust system vented to the outside of the building without recirculation. Each BSC shall have one duct and one exhaust fan. Class II Type B1 BSCs should be installed if routine use of

minute quantities of volatile toxic chemicals, volatile radionuclides, or noxious specimens or cultures is expected.

- Include a leak-tight duct, a damper in the duct near the cabinet to permit flow, adjustment closure and decontamination, and an external exhaust fan as the final system component. **A butterfly valve or air control roll valve (depending on ceiling height) must be installed to control airflow.** An alarm must be installed to indicate loss of exhaust flow. Each B1 BSC should have its own dedicated exhaust system.
- The exhaust system shall meet the exhaust requirements as specified by the BSC manufacturer. Building duct static pressure must be maintained at 1.5 inches water gage (w.g.) (38.1 mm w.g.).
- Provide an alternate room exhaust duct for use when the BSC is "off". A bypass damper using the BSC exhaust fan is acceptable.
- **A Bag-In/Bag-Out Unit for both the exhaust and recirculation is required in B1 BSCs.** This protects the technician handling the filter as well as the environment. See Drawing – Bag-In/Bag-Out Unit.

Type B1 Laminar Airflow BSC Specifications

- The BSC must be National Sanitation Foundation (NSF Standard 49) certified for Class II Biological Safety Cabinets and bear the NSF seal.
- The BSC must be a bench top model.
- Installation of **Ultraviolet (UV) lamps** is not permitted in BSCs at the NIH. The NIH, CDC, NSF/ANSI, and the American Biological Safety Association all agree that ultraviolet (UV) lamps are not recommended nor are they necessary. NSF Standard 49, the industry testing standard for all biohazard cabinetry, does not provide any performance criteria for UV lighting and specifically states in Section 4.24.2 "UV lighting is not recommended in Class II (laminar flow) biohazard cabinetry." Numerous factors affect the activity of the germicidal effect of UV light, which require regular cleaning, maintenance and monitoring to ensure germicidal activity. In addition, there are safety hazards associated with UV light exposure, which include cornea burns and skin cancer.
- Maintain Minimum airflow velocity of 100 ft/ min (0.5 m/s) through the work access opening.
- All biologically contaminated plenums are negative to the room or surrounded by negative pressure plenums.
- Two hose cocks on each side wall (total of four hose cocks) installed and labeled for vacuum in the four and six foot models. No gas cocks can be installed. Open flames disrupt the laminar flow in BSCs; therefore the NIH does not recommend the use of gas. If it is determined by the DOHS that gas must be used in a BSC,

- Electrical supply should be 115 or 230 volt, 50 or 60 Hertz, AC, three prong convenience duplex receptacle for grounded plugs installed in each side wall in the four and six foot models. It should be protected with thermal circuit breakers.
- Five (5) amp breakers shall be utilized on the B1 BSC.
- The BSC must be listed with Underwriters Laboratories to meet requirements of UL 3101-1.
- The BSC must have a self contained electronic module control panel.
- Hinged or sliding front window with audible and visual alarms, which indicate, window is not at the correct position for sliding window models.
- A front access opening of eight (8) inches or ten (10) inches.
- Standard size HEPA filters with keyed, sealed gaskets are readily available from any of the major filter manufacturers.
- HEPA filtered downflow air composed of uncontaminated recirculated inflow air.
- The BSC must be designed so service can be performed (including decontamination) without disconnecting utility services and/or moving the cabinet.
- The BSC shall be ergonomically designed for maximum user comfort and adjustability.
- Noise levels are not to exceed 67 dbA with a maximum background level of 57dbA.
- As shipped, can pass through a standard doorway (6' 7"high x 2' 11" wide) or removable top mounted parts are removed before shipping and labeled with the NIH purchase order number and the manufacturer's serial number, even if parts are shipped on the same pallet.

Two copies of instruction/use booklets furnished with each cabinet and listing:

- Instructions for installing, operating and performing preventative maintenance on the cabinet.
- Identification of service parts with manufacturer's part number, and quantify required for preventative maintenance purposes.
- The fan performance curve for the integral fan(s).
- Results of the required NSF tests for each cabinet purchased.
- Results of all the production quality control tests on each cabinet.


A standard stand made by the same cabinet manufacturer which can be securely fastened to the cabinet shall be ordered for each cabinet, unless otherwise requested by

the end user. The cabinet and stand combination (referred to as a console by NSF) must have passed the NSF “tip test” section A.8.3.1 of *NSF/ANSI 49 2002* when mounted on the stand to be used with the cabinet, in order to be acceptable to the NIH.

Contact the DOSH at 301 496-3457 in Building 13, Room 3K04 or Consult NSF International Standard/American National Standard 49-2002 if you have further questions.

FIGURE 1: The Class II, Type B1 BSC (classic design).


Connection to building exhaust system required.

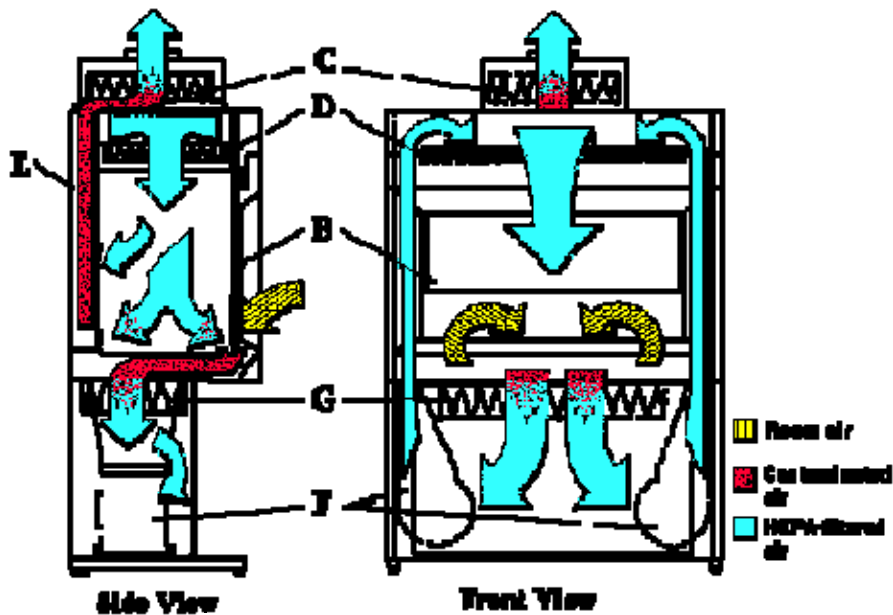

- A. front opening
- B. sash
- C. exhaust HEPA filter
- D. supply HEPA filter
- E. negative pressure exhaust plenum
- F. blower
- G. additional HEPA filter for air supply

Note: The cabinet exhaust needs to be connected to the building exhaust.

FIGURE 2: BAG-IN / BAG-OUT FILTER

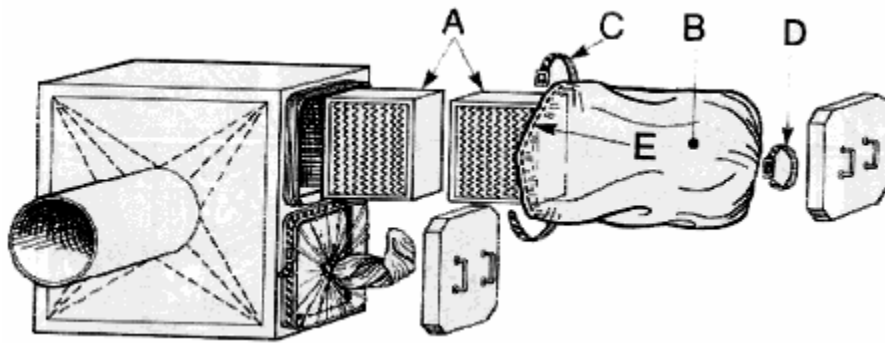

- A.** filters
- B.** bags
- C.** safety straps
- D.** cinching straps
- E.** shock cord located in the mouth of the PVC bag restricts the bag around the second rib of the housing